




GPR108 Antibody

CATALOG NUMBER: 48-150

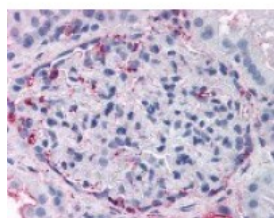

Immunohistochemistry staining of
GPR108 in kidney tissue using GPR108
Antibody.

Specifications

SPECIES REACTIVITY:	Gorilla, Human
TESTED APPLICATIONS:	IHC
APPLICATIONS:	GPR108 antibody can be used in ELISA, Western Blot, and immunohistochemistry starting at 5 ug/mL.
USER NOTE:	Optimal dilutions for each application to be determined by the researcher.
SPECIFICITY:	BLAST analysis of the peptide immunogen showed no homology with other human proteins.
IMMUNOGEN:	GPR108 antibody was raised against a peptide located in the N-Terminal extracellular domain of GPR108 (Human).
HOST SPECIES:	Rabbit

Properties

PURIFICATION:	Immunoaffinity Chromatography
PHYSICAL STATE:	Liquid
BUFFER:	PBS, 0.1% sodium azide.
STORAGE CONDITIONS:	GPR108 antibody should be stored long term (months) at -80 °C and short term (days) at 4 °C. As with all antibodies avoid freeze/thaw cycles.
CLONALITY:	Polyclonal
CONJUGATE:	Unconjugated

Additional Info

ALTERNATE NAMES:	GPR108, G protein-coupled receptor 108, LUSTR2, Protein GPR108
ACCESSION NO.:	Q9NPR9
PROTEIN GI NO.:	296439338
OFFICIAL SYMBOL:	GPR108
GENE ID:	56927

Background

BACKGROUND: LUSTER2 is an Orphan-U GPCR with an unknown ligand.

FOR RESEARCH USE ONLY

December 13, 2016